
Nosebud⁺

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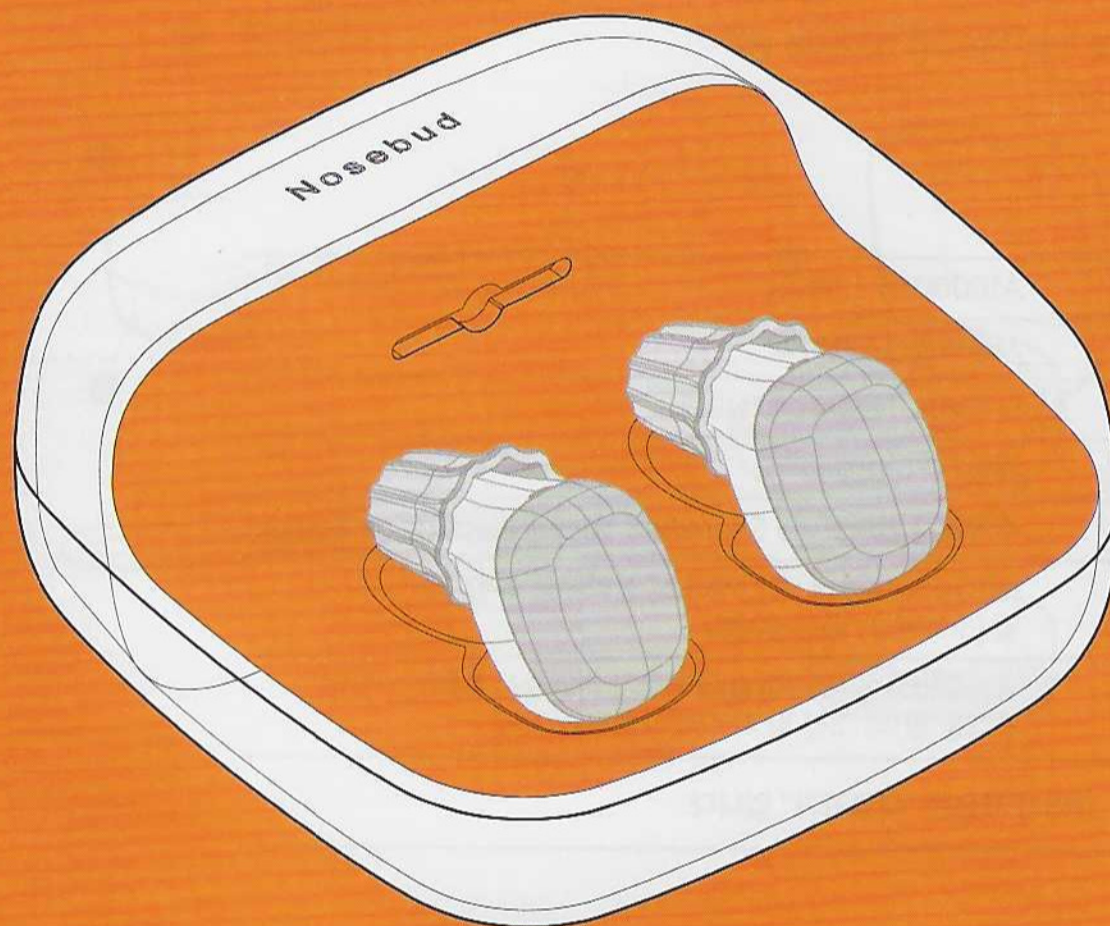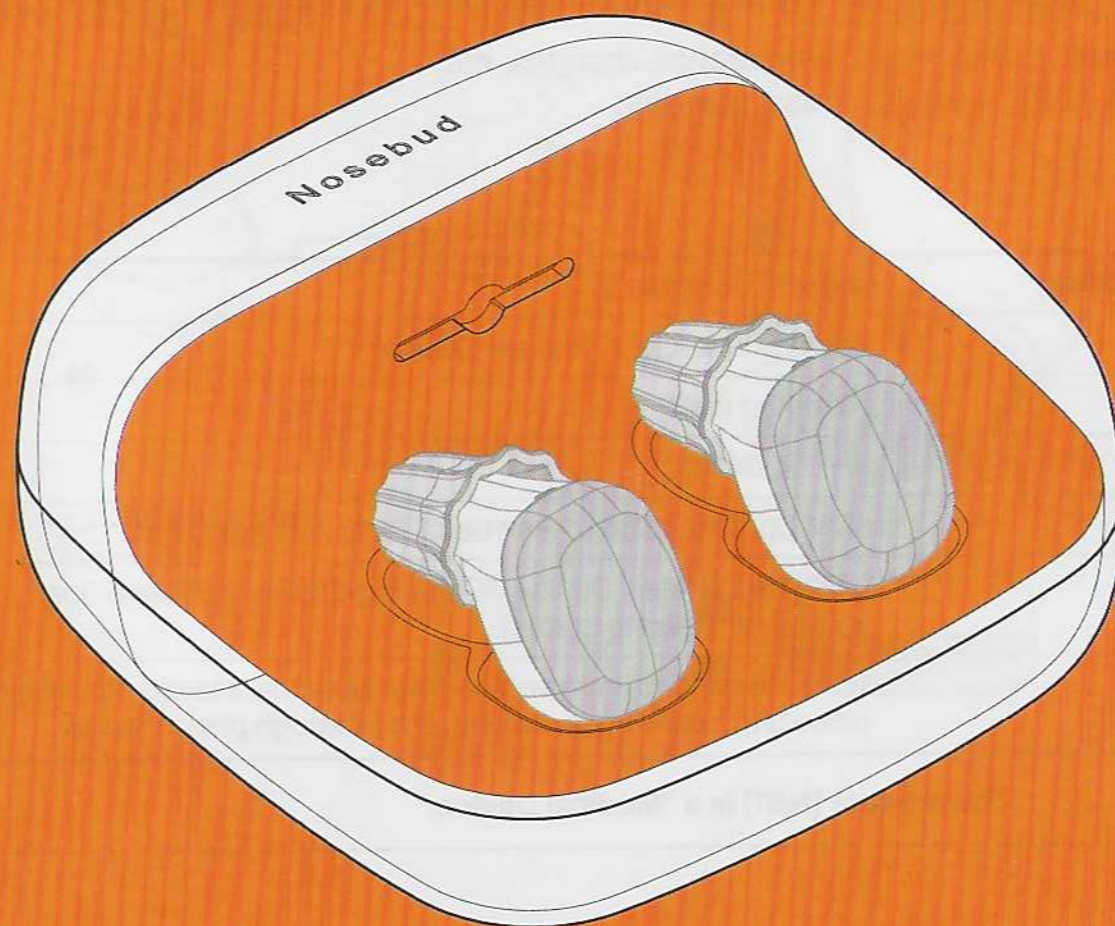
User Manual / Warranty

01 Efficacy of light therapy device

Relieves seasonal allergic rhinitis symptoms by irradiating 660nm red wavelength and 940nm near-infrared wavelength to the rhinitis-affected area in the nose.

| 660nm red wavelength : Antioxidant effect, relief of rhinitis symptoms

The red LED wavelength of 660 nm suppresses free radicals and strengthens the vitality of mucosal dermal cells, alleviating rhinitis symptoms.

| 940nm near-infrared wavelength : Suppress allergic reactions, relieves rhinitis symptoms

The near-infrared wavelength of 940 nm suppresses allergic reactions by inducing appropriate histamine production by mast cells.

| Allergic rhinitis : Hypersensitivity of the nasal mucosa to certain substances

Allergic rhinitis is a disease in which inflammatory cells flock to the irritated area after the nasal mucosa is exposed to allergen-causing substances, causing an inflammatory response due to various mediators secreted by them.

- Causes: pet hair, fine dust, pollen, mites, temperature changes, etc.
- Symptoms: runny nose, stuffy nose, itchiness, sneezing

02 Product overview

Medical device certification [class]	A16060.02 [class 2]
Class name	Lowpowered Light Therapy Device
Product name [Model]	Nosebud [NB1]
Certification number	제인24-648호
manufacturer	NTV Co., Ltd. / Shenzhen Yuefull Electronics Co., Ltd
Packaging unit	[NB1]1Set (2 Main Device, 1Charging Cradle, 4 Silicone tips) and 1 [UV charging cradle]
Purpose of use	A device used to improve and treat seasonal allergic rhinitis symptoms (sneezing, runny nose, itching, etc.) through phototherapy.

The product [NB1] is a “Medical Device.”

03 Specifications

Internal power	DC 3.7V
Wavelength (frequency)	660nm & 940nm
Output (style)	Continuous operation
Operating time	Around 3 minutes
Safety function	1. Automatically turns off after 3 minutes of operation. 2. Touch the power switch for 3 seconds to turn off the power.
Electric shock protection	Internally powered devices, BF type devices


04 Product composition

[NB1] 2 Nosebud main device, 1 Portable charging cradle, 1 Silicone tip container, 4 Silicone tips (2 regular, 2 slim) in the container, 1 Instruction manual


[UV Charging cradle] 1 Charging Cradle (Charging cable is not included.)

A. Nosebud [NB1]

- 1 Main device
- 2 Light guide
- 3 Power switch
- 4 Indicating light


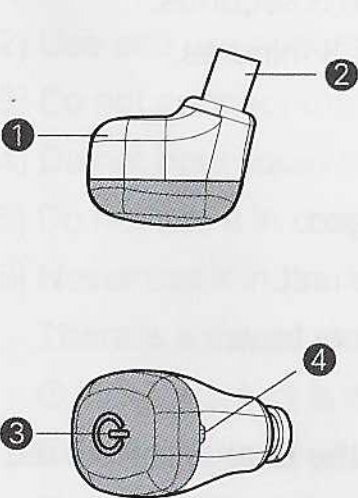
- 5 Portable charging cradle
- 6 Indicating light
- 7 USB-C charging terminal
- 8 Lanyard

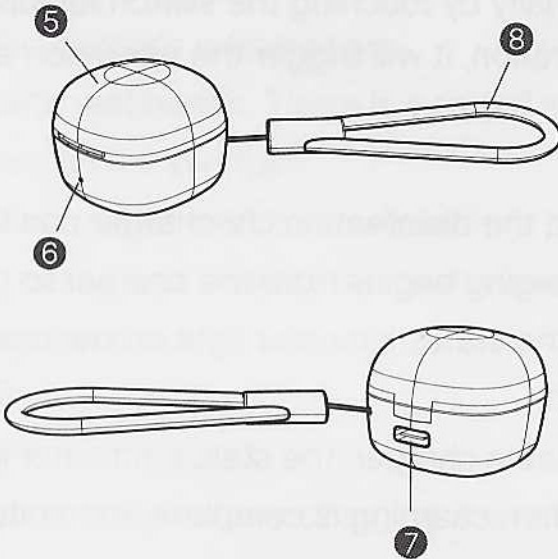B. UV charging cradle

- 9 UV charging cradle
- 10 UVC sterilizing lamp
- 11 Indicating lights
- 12 USB-C charging terminal

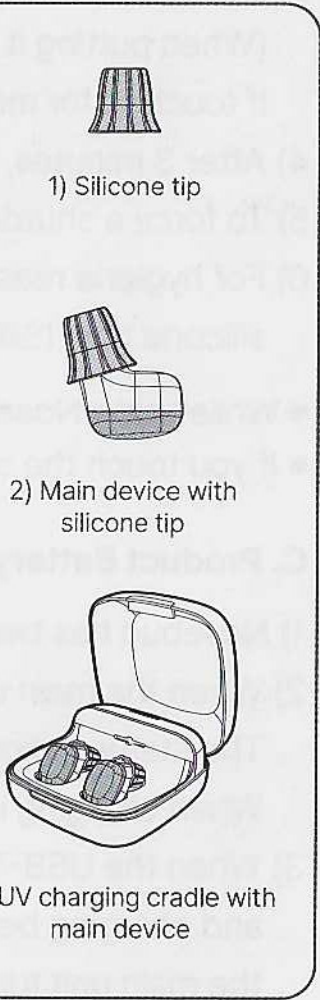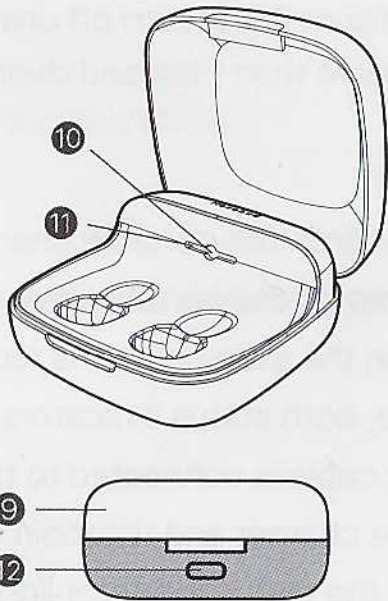05 How to use


A. Preparation before use

- 1) Before turning on the power, check the location and environment are appropriate.
- 2) Do not leave liquids around the device that may spill and affect the device.
- 3) Since flammable liquids and gases may ignite at the electrodes due to electrical sparks when charging, be sure to check for flammable gases before using the device and be careful.
- 4) Avoid areas with strong electricity, magnetic fields, or areas that may be exposed to heat or humidity.
- 5) Check if the main device (low-power light irradiator) is operating normally.
- 6) Never use this product for purposes other than its intended use.

B. How to use and operating sequence

- 1) Place the enclosed silicone tip to the light guide. (There are two types: regular and slim.)
- 2) Touch the main device power switch for 1 second, the light guide turns on for 3 minutes.
- 3) Gently push the silicone tip side into your nasal cavity.
(When putting it in your nose, make sure not to touch the switch for more than 3 seconds. If touched for more than 3 seconds, the power may be turned off.)
- 4) After 3 minutes, the light guide automatically turns off.
- 5) To force a shutdown, touch the power switch for 3 seconds.
- 6) For hygiene reasons, please do not share silicone tips with others and use separate silicone tips. (Silicone tips can be purchased additionally.)

- ※ While using Nosebud, make sure not to turn off unintentionally by touching the switch for longer than 3 seconds.
- ※ If you touch the switch for more than 1 second during operation, it will trigger the operation again for 3 minutes.

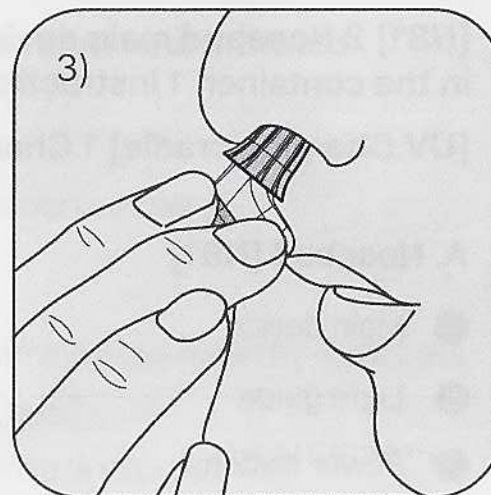C. Product Battery Charging

- 1) Nosebud has two chargers. Both the portable charger and the disinfection UV charger can be charged.
- 2) When the main unit is placed on the portable charger, charging begins from the charger to the main unit.
The status indicator light on the main unit turns red, and the status indicator light on the charger turns blue.
When charging is complete, both status indicators turn off.
- 3) When the USB-C charging cable is connected to the portable charger, the status indicator light on the charger blinks red and charging begins for the charger and the main unit. When charging is complete, the status indicator light on the main unit turns off and the status indicator light on the charger stops blinking and turns on.
- 4) When the main unit is placed on the disinfection UV cradle, charging begins with the main unit once the lid is closed.
The status indicator light on the main unit turns red. When charging is complete, the status indicator light turns off.
- 5) When the USB-C charging cable is connected to the disinfection stand charger, the status indicator light on the charger blinks red and charging begins for the charger and the main unit. When charging is complete, the status indicator light on the main unit turns off and the status indicator light on the charger stops blinking and turns on.
- 6) The charger for the disinfection stand has a built-in UVC ultraviolet light. When you close the charger lid, the blue light turns on and the UVC disinfection light starts working. The UVC disinfection light and blue light turn off when you open the lid and turn on again when you close it. The UVC disinfection light and blue light automatically turn off after 3 minutes.
(Do not look directly at the UVC disinfection light.)

D. Storage and management method after use

- 1) To reuse the used silicone tip, remove it from the product and wipe it with an alcohol swab or wash it in water and dry it completely before using.

- 2) When cleaning the device after use, be sure to disconnect the charging cable from the device and wipe off any dust, sebum, etc. on the surface of the device with a soft, dry cloth or tissue.
- 3) The device may become deformed, so store it at room temperature and avoid high temperature and humidity.
- 4) Do not store it in a place where chemicals are stored or gas is generated. 5) When not in use for a long period of time, be sure to separate the product and charger.
- 6) Store in a safe place away from vibration or impact.
- 7) Keep out of reach of infants.

06 Precautions for use

A. Warning

- 1) Do not disassemble/repair/modify the product arbitrarily.
- 2) Use and use according to the doctor's prescription.
- 3) Do not connect the power with wet hands. There is a risk of electric shock and fire.
- 4) Do not pour liquid substances on the product.
- 5) Do not use it in conjunction with other medical devices.
- 6) Never use it in the following cases and contact the customer service center for inspection and service.
There is a risk of electric shock and fire.
 - ① If the product is damaged
 - ② If the product gets wet
 - ③ If the USB connector and cable for battery charging are damaged
- 7) The ME device must be connected to an appropriate power source.
(When using an adapter, check the KC certification and rating (Input: 100-240VAC, 50/60Hz, 0.35A, Output: 5Vdc, 1A) before use.
- 8) If the ME device is not to be used for a certain period of time, remove the battery. However, battery removal is performed only by our company and the service provider designated by our company.
- 9) The light emitter is a light guide and silicone tip. Please check when using the product.

B. Precautions regarding the age, gender, or health of the target user

- 1) All ages can use it. However, if your skin is sensitive to light, do a small sensitivity test on your sensitive skin before use, and patients with photosensitivity should not use this product. Also, pregnant women, lactating women, women of childbearing age, newborns, infants, children, and the elderly are prohibited from using it.
- 2) Do not allow children to touch or use the device alone.
- 3) If possible, the target of phototherapy should undergo an endoscopic examination to rule out polyposis and sinusitis.
- 4) The following people should consult a doctor before using it.

- Patients with acute diseases
- People with abnormal skin in the area of use
- Patients with high fevers
- Patients with malignant tumors
- People with sensory impairment

C. Precautions regarding adverse reactions to medical devices that may occur as a result of using the medical device, fatal side effects, and accidents due to carelessness in use

- 1) Infrared rays are invisible to the eye. Do not put your eyes close to the infrared LEDs to check if they are working.
- 2) In the following cases, be sure to consult a doctor before using the product.
 - ① If you have or are suspected of having any kind of disease (patients with malignant tumors/cancer patients/patients with high fever/people with extremely weakened vitality/users of life-extending devices, etc.)
 - ② If you have an allergy or hypersensitivity to heat when touching your skin
 - ③ Patients who have recently undergone surgery, pregnant women.
 - ④ If your skin that the product touches has inflammation or trauma, etc.

D. General precautions

- 1) Do not forcefully wrap or bend the cable. There is a risk of malfunction, electric shock, or fire.
- 2) Do not place flammable substances such as candles, cigarettes, or alcohol, or containers containing water, heavy objects, or sharp objects on top of the product. This may cause electric shock, fire, or damage.
- 3) Do not impact or drop the product.
- 4) Do not clean the device with organic solvents (benzene, thinner, etc.).
- 5) Keep the product away from heating appliances such as stoves.
- 6) Do not use excessive force on the power button. It may cause product failure.
- 7) Do not hold the charging cable and swing or turn it while it is connected. It may cause injury and product failure.

E. Precautions for interaction

- 1) Be careful as there is a risk of danger when using it with other devices
- 2) Do not use it if you have an allergic reaction to light or are taking medications that may cause such reactions.
- 3) Do not use it with objects that may cause electric shock.
- 4) Do not use it with alcohol or flammable liquids/gases.

F. Precautions for use by pregnant women, nursing women, women of childbearing age, newborns, infants, children, and the elderly

- 1) Pregnant women, nursing women, women of childbearing age, newborns, infants, children, and the elderly should refrain from using it as their skin may be sensitive and may cause allergies.

G. Caution for application

- 1) Do not expose the product to excessive heat such as sunlight or fire.
- 2) Do not apply skincare products, cosmetics, or medicines to the nose and around the nose before and while using the product.
- 3) If you have sensitive skin, consult a doctor before using.
- 4) If you feel hot, immediately turn off the product by touching the power button to avoid the risk of burns.
- 5) Do not bend or fold the device severely by hand. This may cause damage or breakage of the product.
- 6) When disposing of this product, please dispose of it safely in accordance with the recycling regulations of your area.
Do not disassemble the product when disposing of it.

H. Precautions for preventing safety accidents

- 1) Keep away from heat.
- 2) Do not place or use the product near flammable or combustible materials (insecticides, hairspray, gas, oil, or other easily flammable items) or in humid places.
- 3) Check the device frequently to see if it is worn out, broken, or damaged.
- 4) Turn off the power when not using the device.
- 5) Do not use it in places with high humidity such as bathrooms, and use it indoors.
- 6) Do not subject the product to strong impacts or place heavy objects on top of the device. 7) This product must be used by the patient himself/herself.

07 Storage method

A. Storage conditions or storage method

- 1) Store in a place that does not come into contact with water.
- 2) Store in a place free of dust and moisture.
- 3) Store with care to maintain stability, such as tilt, vibration, and impact (including during transport).
- 4) Do not store in a place where chemicals are stored or gas is generated.
- 5) Store in a place away from direct sunlight.

B. Environmental conditions

- 1) Operation: Temperature (0°C~35°C), humidity (20%~80%), atmospheric pressure (800~1060hPa)
- 2) Storage: Temperature (-10°C~40°C), humidity (20%~80%), atmospheric pressure (800~1060hPa)

Product Warranty

Product Name		Nosebud	
Model		NB1	
Date of purchase			
Warranty		1 year after purchase	

customer	name		Phone no.
	address		
store	name		Phone no.
	address		

Please be sure to fill in the blank information.

Consumer Compensation Information

Consumer Damage Type			Compensation details	
			Within the warranty	After the warranty period
In case of performance or functional failure that occurs under normal use conditions (Parts retention period within 3 years)	If major repairs are required within 10 days of purchase		Product exchange or refund	
	If repair of important parts is required within 1 month of purchase		Product Exchange	
	If the exchanged product requires major repairs within 1 month		refund	
	When exchange is not possible			
	repair possible	In case of defect	Free repair	Paid repair
		The integrated defect was repaired, but the problem reoccurred		Product exchange or refund
	If repair is not possible			
	If repair is not possible because repair parts are not available			
	If the business owner loses the product that the consumer requested repairs			
Damages incurred during product purchase, transportation, product installation, etc.		Product exchange		
functional failure due to consumer	When repair is possible		Free repair	Paid repair
	When repair is not possible		Product exchange after collecting paid repairs	company's standards
a malfunction occurs due to external causes other than a defect in the product itself or unauthorized repair			Paid repair	Paid repair
The above compensation provisions are in accordance with Article 12, Paragraph 2 of the Consumer Protection Act. Matters not specified are subject to consumer damage compensation regulations. This product is covered by production and sales liability insurance.			Customer Service Center : Please refer to the seller's website	